



PCT/GB 2003 / 000974

10/506906 #2
Rec'd PCT/PTO 17 SEP 2004

INVESTOR IN PEOPLE

**PRIORITY
DOCUMENT**SUBMITTED OR TRANSMITTED IN
COMPLIANCE WITH RULE 17.1(a) OR (b)The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

REC'D 02 MAY 2003

WIPO PCT

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.


In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed

Dated 1 April 2003

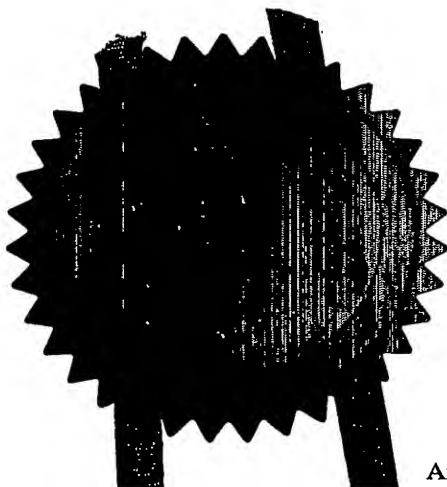**BEST AVAILABLE COPY**

An Executive Agency of the Department of Trade and Industry

07-03-02

04:07pm

From-D YOUNG & CO

+02380 719800

T-511 P.002


F-513

Patents Form 1/77

Patents Act 1977
(Rule 16)07MAR02 E701648-1 002246
P01/7700 0.00-0205394.0

Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

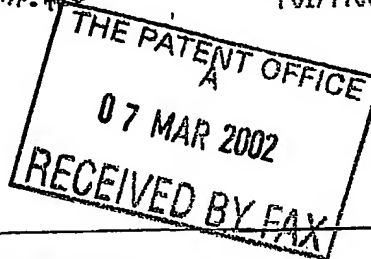

The Patent Office

Cardiff Road
Newport
South Wales
NP10 8QQ

1. Your reference

P003500GB ATM

2. Patent application number

(The Patent Office will fill in this part)

0205394.0

07 MAR 2002

3. Full name, address and postcode of the or of each applicant (underline all surnames)Medical Research Council
20 Park Crescent
London W1B 4AL
United Kingdom

"SEE CONTINUATION SHEET"

Patents ADP number (*if you know it*)

5840624001

If the applicant is a corporate body, give the country/state of its incorporation

United Kingdom

4. Title of the invention

sCD Molecule

5. Name of your agent (*if you have one*)

D Young & Co

"Address for service" in the United Kingdom to which all correspondence should be sent (*including the postcode*)21 New Fetter Lane
London
EC4A 1DAPatents ADP number (*if you know it*)

59006 ✓

6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (*if you know it*) the or each application number

Country

Priority application number
(*if you know it*)Date of filing
(*day / month / year*)

7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application

Number of earlier application

Date of filing
(*day / month / year*)8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (*Answer 'Yes' if:*

Yes

- a) any applicant named in part 3 is not an inventor, or
- b) there is an inventor who is not named as an applicant, or
- c) any named applicant is a corporate body.

See note (d))

Patents Form 1/77

0030819 07-Mar-02 04:44

Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form. Do not count copies of the same document

Continuation sheets of this form

Description 33 ✓

Claim(s) 3 ✓

Abstract 1 ✓

Drawing(s) 8 only ✓

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (Patents Form 1/77)

Request for preliminary examination and search (Patents Form 9/77)

Request for substantive examination (Patents Form 10/77)

Any other documents Appendix 1
(please specify)

11.

I/We request the grant of a patent on the basis of this application.

Signature

D Young & Co

Date 07 March 2002

D Young & Co (Agents for the Applicants)

12. Name and daytime telephone number of person to contact in the United Kingdom

Antonio Maschio

023 8071 9500

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 08459 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

Patents Form 1/77

CONTINUATION SHEET**OUR REF: P008500GB ATM****3. Further Applicants**

Addenbrookes NHS Trust
Hills Road
Cambridge
CB2 2QQ

8 339 632 001

County/state of incorporation United Kingdom

Cambridge University Technical Services Limited
The Old Schools
Trinity Lane
Cambridge
CB2 1TS

6956809004

County/state of incorporation United Kingdom

sCD MOLECULE

The present invention relates to the use of cluster of differentiation (CD) molecules in detecting the presence and/or assessing the progression of one or more disease states in an individual. In particular it relates to the use of profiles/fingerprint/s of secretory CD (sCD) molecules in detecting and/or assessing the progression of one or more disease states in an individual. Further uses of sCD profiles according to the present invention are also described.

Background to the invention

10

Rapid and accurate diagnosis is essential in medicine as in many cases early diagnosis and successful treatment correlates with a better outcome and reduced hospitalisation. Currently, the clinical diagnosis and staging of many diseases of global significance involve different invasive procedures such as histopathological analysis of biopsy samples which are usually obtained when the disease process is at a relatively advanced stage. In many cases, a classic histopathological approach may not be sufficient to produce accurate diagnosis and any delay in confirming the diagnosis would have financial and morbidity repercussions for the healthcare institution and most importantly for the individual. Disease states and disease staging are also determined by different imaging techniques such as X-rays, nuclear magnetic resonance (NMR), CT analysis and others, however, these are expensive and impractical when dealing with large numbers of samples.

A variety of diseases have a genetic basis and can be characterised by changes in the overall patterns and/or expression levels of various proteins. For example, some cancers are associated with changes in the expression of oncogenes or tumour suppressor genes. Furthermore, disease conditions or disorders associated with dysregulated cell cycle and development can be attributed to changes in transcriptional regulation of specific genes.

30

Although there are several genetic assays available to assess gene mutations, the identification of specific genetic changes may not always be a direct indicator of a disease or disorder and thus can not be relied upon as an accurate prognostic indicator.

- 5 Certain genetic changes are exhibited by alterations in cell surface antigens. Again, however, prior attempts to develop a diagnostic assay for complex disease conditions or disorders based on the identification of single antigen have not been uniformly successful.
- 10 In addition, or alternatively, biochemical analysis of a patient may be used to diagnose disease state. For example, the presence of Bence Jones protein in urine is an indicator that an individual has multiple myeloma. However, classical biochemical methods are limited for example an elevated cholesterol in serum indicates hypercholesterolaemia but does not, definitively indicate atherosclerosis. A further disadvantage of
- 15 biochemical methods of diagnosis is that they generally permit the measurement of only one or two indicator/s of disease in any one test. Consequently, they provide an incomplete picture of the disease state of an individual. Moreover, if several tests are performed in an attempt to provide a more complete picture, this inevitably increases the number of variables which complicates interpretation. Furthermore, for many
- 20 diseases there are no reliable biochemical markers, especially for diseases of global importance such as breast cancer, colorectal cancer and lung cancer. In the case of colorectal cancer, a number of carcinoembryonic antigen (CEA) markers have been identified, however they have poor sensitivity and very low specificity. The situation is similar with disease conditions requiring surgical intervention. There is still, for
- 25 example, no marker for acute appendicitis and consequently, a great many patients undergo unnecessary invasive surgery. It has been estimated that more than 40,000 unnecessary appendicitis operations occur each year due to misdiagnosis with associated costs of \$700 million. In a recent larger retrospective study, Flume and colleagues show that misdiagnosis occurs in 15% of instances.

Therefore, there is a pressing need in the art to provide a simple and complete picture of the disease state or condition of an individual. Such a 'picture' would be of use in predicting and/or detecting the presence of a disease or condition, in assessing the therapeutic strategies and the potential of various agents and in monitoring the progression and successful treatment of disease states or condition.

Lymphocytes and other leukocytes express a large number of different antigens associated with their outer plasma membranes that can be used to identify distinct functional cell subsets. Many of these antigens were "classically" known to be receptors for growth factors, cell-cell interactions and immunoglobulins; molecules for cell adhesion or complement stimulation; enzymes and ion channels. A single systematic nomenclature has been adopted to classify monoclonal antibodies to human leukocyte cell surface antigens termed cluster of differentiation (CD) antigens, also referred to as 'CD molecules (Kishimoto *et al.*, 1996 Proceedings of the Sixth International Workshop and Conference held in Kobe, Japan. 10-14 Garland Publishing Inc, NY, USA).

The data required in order to define a CD has changed over the years, not surprisingly in view of the advances in modern technology. Initially, clustering depended absolutely on the statistical revelation of similarities in reaction pattern of two or more antibodies, analysed on multiple tissues. It is now accepted that CD molecules are to be classified by molecular characteristics. Thus it has become customary to use the CD marker (for example CD21) to indicate the molecule recognised by each group of monoclonal antibodies. The current list of CD markers is constantly updated as new antigens are identified and eventually, the CD list will encompass all human lymphocyte cell surface antigens and their homologues in other mammalian and non-mammalian species (Mason *et al.*, 2001, *Immunology*, 103, 401-406).

Historically, CD cell surface antigens have been used as markers in diagnosis. Indeed leukemias are diagnosed on the basis of cell morphology, expression of specific CD antigens, lymphoid (LY) and myeloid (MY) antigens, enzyme activities and cytogenetic abnormalities such as chromosome translocations. The expression of up to

three CD antigens on leukemia cells is determined using labelled antibodies to particular CD antigens with analysis by flow cytometry.

5 Significantly, however, it has been observed that often (if not always in normal or disease states) the surface bound CD immunological specificity molecules (intact CD molecules or fragments thereof) are found soluble in the serum and in other body fluids. Subsequent research has shown that indeed CD molecules can be secreted from cells as a result of "active" processes such as alternative splicing (Wolfson and
10 Milstein, PNAS, 91 (14) 6683-6687) or "passive" processes such as cell surface shedding. Thus, CD molecules can be found in two forms, membrane associated CD molecules and secretory CD molecules (sCD) which can be complete molecules or fragments thereof.

15 It is generally accepted however that the change in levels of any one sCD is not indicative of any given disease state and cannot therefore be used in diagnosis of disease states.

Recent studies (those of WO 00/39580) have described a system for the diagnosis of leukemia, whereby immunoglobulins are immobilised on a solid support and are used
20 to detect cell-surface antigen levels, in particular cell-surface CD antigen levels in samples of cells. Using this approach, a pattern of expression of cell-surface bound CD antigens is generated which the inventors have shown to be indicative of the presence of various defined leukemias in a patient. However, there are several disadvantages with this technique. Firstly and importantly, it is a cell-based technique. Such
25 techniques have many disadvantages associated with them, for example that of background noise and the difficulty of measuring antigen levels accurately. In addition, generally, it is only practical to test one cell type at a time (due to the problems of background noise and providing a suitable non-diseased control).

30 Therefore, there still exists a need in the art for a simple method for diagnosis of different diseases and conditions by the measurement of CD antigens wherein such method produces a complete and accurate picture of disease.

Summary of the invention

The present inventors have surprisingly found that particular disease states can be characterised by specific patterns of levels of secretory CD molecules derived from the body fluids of an individual. That is, the profile of the levels of sCD antigens correlate with particular diseases or disorders. This finding is especially surprising since the levels of sCDs found in the body fluids of an individual are generally very low, and as the sCD released by cells would only be expected to change in some, and not all cell types of an individual when affected by one or more diseases, then the change of levels of secreted CD levels detectable in the body fluids of diseased individuals as compared with non-diseased individuals would be expected to be minimal.

Thus, in a first aspect, the present invention provides a secretory CD (sCD) fingerprint of one or more disease states.

In the context of the present invention, the term 'CD' refers to a different cell surface leukocyte molecule recognised by a given monoclonal or group of monoclonal antibodies which specifically 'cluster' to the antigen/molecule in question. Many, if not all of these molecules produce secretory forms by alternative splicing, shedding or other mechanisms. Thus in the context of the present invention, the term 'secretory CD molecule (sCD)' refers to a secreted form of a cell surface leukocyte molecule recognised by a given monoclonal or group of monoclonal antibodies as herein described.

As herein defined, the term 'secretory CD fingerprint (sCD)' describes the pattern or profile of levels of more than one secretory CD in one or more individuals. A sCD fingerprint as herein defined may be representative of one or more non-diseased individual/s or a one or more diseased individual/s. Preferably, a secretory CD fingerprint describes the level of five or more secretory CD molecules, more preferably it is 6, 7, 8, 9, 10 or more sCD molecules, more preferably still a secretory

CD fingerprint describes the levels of 15 or more sCD molecules. Most preferably a secretory CD fingerprint describes the levels of 20 or more sCD molecules.

- 5 sCD levels from normal and/or diseased individuals may be collated in order to generate one or more reference sCD fingerprints. A 'reference' sCD fingerprint is a fingerprint which is advantageously generated from sCD measurements from more than one individual and is representative of the sCD levels of either a 'normal' or diseased individual. Advantageously, these reference fingerprints are collected together to form a database so that abnormal fingerprints generated from patient
- 10 samples can be distinguished from normal reference fingerprints in the database. In addition, by comparing one or more patient sample fingerprints with one or more reference fingerprint/s corresponding to one or more diseases, then the disease state of an individual may be established.
- 15 A secretory CD (sCD) fingerprint may be generated from one individual. Preferably however, each fingerprint is generated from more than one individual. Advantageously, it is generated from more than five, ten, fifteen or twenty individuals. One skilled in the art will appreciate that the greater the number of individuals used to generate the reference fingerprint, then the more representative of any given disease
- 20 state or of a normal individual the reference sCD profile/s will be. Fingerprints may be simplified by using the average values of the data obtained for each sCD for a number of individuals. For example, the modal value is used for data obtained from a very small number of samples.
- 25 One skilled will appreciate that often more than one disease state may be present in an individual at a given time. This may complicate the CD fingerprint obtained, such that the fingerprint is an aggregate fingerprint of several disease states. The effect of multiple disease states (composites) in an individual may be minimised if the reference fingerprint for any given disease state is generated from several or many individuals.
- 30 Importantly, composites may generate their own patterns and be used as reference in their own right.

Measuring the levels of sCD molecules is carried out using methods known to one skilled in the art and described herein.

5 sCD levels are measured in samples of body fluids. Suitable body fluids for measuring sCDs include any one or more selected from the group consisting of the following: tissue fluid, serum, blood, cerebrospinal fluid, urine, lymphatic fluid, aspirate, bone marrow aspirate and mucus. One skilled in the art will appreciate that this list is not intended to be exhaustive.

10 In a further aspect the present invention provides a method of generating a secretory CD (sCD) fingerprint of one or more disease state/s comprising the step of measuring the levels of more than one secretory CD in one or more samples from one or more individuals and collating the data.

15 According to the above aspect of the invention, the samples for testing and used to generate a given sCD profile may be from one body fluid type or more than one body fluid type. Advantageously the one or more samples for testing and used to generate a sCD profile are taken from more than one body fluid for any given disease state.

20 As used herein, the term 'collating' the data means to put the data into a form so that one or more pattern/s of the levels of sCDs within that disease state is apparent.

25 As the present inventors surprisingly found that a sCD fingerprint of a diseased individual is different from that of a non-diseased individual, it was realised that by comparing the sCD fingerprint/profile of a sample from a diseased patient with that of one or more reference sCD fingerprints representing one or more defined disease states then the presence and nature of a disease in that individual could be ascertained.

30 Thus, in a further aspect, the present invention provides a method for predicting the presence of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or more reference sCD fingerprint/s.

In the context of the present invention, the term 'predicting the presence of one or more disease states' refers to the process of detecting the presence of one or more disease states before the onset of the clinical signs of the disease are apparent in the individual. The clinical signs of disease are characteristic of each disease state or
5 group of disease states, as long as the disease is present in that individual.

As referred to herein, the comparing step may refer to comparing an individual's sCD fingerprint/profile with one or more reference sCD fingerprint/s of one or more disease
10 state/s and/or with a reference sCD fingerprint of a non-diseased 'normal' individual.

In a further aspect still, the present invention provides a method for detecting the presence of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or
15 more reference sCD fingerprint/s.

The term 'detecting the presence of one or more disease states' refers to the detection of the presence of one or more disease states in an individual once the clinical signs of one or more disease states are apparent in that individual. In addition, the term refers to
20 the process of detecting the presence of one or more disease states in an individual other than the disease state whose clinical signs are apparent in that individual.

According to the above two aspects of the invention, the reference sCD fingerprint/s may be from non-diseased (normal) individuals and/or from diseased individuals.
25 Preferably, they are from diseased individuals.

In a further aspect still, the present invention provides a method for detecting the extent of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or more
30 reference sCD fingerprint/s.

As referred to above, the term 'detecting the extent of one or more disease states' includes within its scope detecting the severity of one or more disease states within an individual. For example it may allow low grade and high grade forms of the disease to be distinguished. It allows localised and metastasised forms of a particular disease to be distinguished. In such cases one or more sCD fingerprints of an individual are compared with one or more reference disease sCD fingerprints representative of disease states at one or more degrees of severity. For example, in the case of neoplastic disease the presence or absence of metastasis may be detected using the method of the present invention.

10

In a further aspect, the present invention provides a method for assessing the progression of a disease state in an individual comprising the step of comparing the sCD fingerprint of an individual at two or more periods during the occurrence of the disease.

15

In the context of the present invention, the term 'assessing the progression of a disease state' means assessing whether the disease has increased in severity, decreased in severity or remains the same severity compared with a different period during the life-span of the disease.

20

The term 'period' in the context of the present invention, generally refers to a time period.

As defined herein, the term a 'disease state' refers to any impairment of the normal physiological functions affecting an organism or any disease condition, disorder or the presence of a particular microbial, viral, parasitic or other pathogenic agent known to one skilled in the art. Suitable disease states for analysis as described herein include but are not limited to: infectious, neoplastic, autoimmune, metabolic, degenerative, psychological, psychiatric, iatrogenic, inflammatory, drug or toxin related, vascular, traumatic and endocrine diseases. Advantageously, 'a disease state' as herein defined refers to any one or more disease selected from the group which includes but is not limited to: infections such as bacterial, fungal, protozoan and viral infections, non-

25

30

neoplastic disorders; stroke; heart condition; atherosclerosis; pain; diabetes, obesity; anorexia; bulimia; asthma; Parkinson's disease; thrombosis; acute heart failure; hypotension; hypertension; urinary retention; metabolic bone diseases such as osteoporosis and osteopetrosis; angina pectoris; hepatitis; myocardial infarction; 5 ulcers; asthma; allergies; rheumatoid arthritis; inflammatory bowel disease; irritable bowel syndrome benign prostatic hypertrophy; pancreatitis; chronic renal failure and psychotic and neurological disorders, including anxiety, schizophrenia, manic depression, delirium, dementia, severe mental retardation and dyskinesias, such as Huntington's disease or Gilles de la Tourette's syndrome and others. Most preferably it 10 refers to appendicitis; Bence Jones' Proteinuria; Chronic Myeloid Leukemia; Colorectal cancer; chronic renal failure; Crohn's Disease; Diabetic Nephropathy; Cardiac pathology; Infection; Liver damage; Lymphoma; macrocytic anaemia; Prostate Cancer; Oligoclonal Banding and Pulmonary Embolism/Deep Vein Thrombosis (eg DVT/PE). One skilled in the art will appreciate that this list is not 15 intended to be exhaustive.

Examples of secretory cluster of differentiation molecules suitable for measurement to generate a sCD fingerprint for use in the methods of the present invention include but are not limited to CD14, CD25, CD31, CD44, CD50, CD54, CD62E, CD62L, CD86, 20 CD95, CD106, CD116, CD124, CD138, CD141, CD40L, CD8, CD23, CD30, CD40. Those skilled in the art will be aware of other suitable sCD molecules for analyses according to the methods of the present invention.

Measuring CD levels may be carried out using methods known to those skilled in the 25 art and described herein. Secretory CD (sCD) levels may suitably be measured in samples of tissue fluids which include, but are not limited to: serum, plasma, lymph fluid, pleural fluid, synovial fluid, follicular fluid, seminal fluid, amniotic fluid, milk, whole blood, urine, cerebrospinal fluid (CSF), ascites, saliva, sputum, tears, perspiration, and mucus. Advantageously, sCD levels are measured from samples of 30 serum using reagents suitable for detecting secretory CDs that include but are not limited to antibodies raised against those CDs. Preferably monoclonal antibodies or engineered antibodies, including phage antibodies raised against secretory CDs or their

membrane bound forms are used for their detection. Commercially available kits for measuring CD levels include those from Diaclone 1, Bd A Fleming BP 1985 F-25020 Besancon Cedex-France and Medsystems Diagnostics GmbH, Rennweg 95b, A-1030 Vienna Austria.

5

Suitable techniques for measuring levels of sCDs include but are not limited to immunoassay including ELISA using commercially available kits such as those described above, flow cytometry particularly multiplexed particle flow cytometry as herein described. Those skilled in the art will be aware of other suitable techniques for measuring CD levels in samples from an individual including antibody 'chip' array type technologies.

10

Secretory CD levels are also measured in a number of individuals with one or more disease states as herein defined, such as appendicitis and the like. Generally one or more secretory CD levels will be elevated in a given disease state as compared with the range found in normal individuals. However, some sCD levels decrease in some disease states compared with the range found in 'normal' individuals and such decreases may also form part of a sCD fingerprint of the present invention. Thus, by measuring the ranges of levels of various secretory CD levels found in a number of individuals with one or more defined disease state/s a 'fingerprint' of secretory CD levels for any defined one or more diseases is generated.

15

20

Preferably, the sCD fingerprint of an individual is generated from any two or more sCDs selected from the group consisting of the following: CD14, CD25, CD31, CD44, CD50, CD54, CD62E, CD62L, CD86, CD95, CD106, CD116, CD124, CD138, CD141, CD40L, CD8, CD23, CD30, CD40. One skilled in the art will appreciate that this list is not intended to be exhaustive and may include CD homologues of human and other mammalian or non-mammalian species. One skilled in the art will appreciate though that in general animal reference sCD fingerprints cannot be used to analyse human diseases and *vice-versa*.

25

30

The invention can also be used for testing human and other mammalian and non-mammalian species using sCD fingerprints from the appropriate animal.

5 One skilled in the art will appreciate that the methods of the present aspect of the invention can be used to test potential therapeutic agents suitable for the prophylaxis and/or treatment of diseases. An agent of therapeutic potential will affect the sCD profile or 'fingerprint' of the disease: If several fingerprints are taken at various stages of a disease and compared with those obtained from samples in which an individual has been treated with a potential therapeutic agent, then the effect on one or more sCD
10 fingerprints can readily be assessed.

In addition, the method of the present invention may also be used to monitor patient compliance with taking a particular drug (agent), and/or undergoing a particular treatment regime.

15

Thus, in a further aspect still, the present invention provides a method for assessing the effect of one or more agent/s on one or more disease states in an individual comprising the step of comparing a sCD fingerprint of an individual at two or more different time periods.

20

According to the above aspect of the invention, preferably the agent is a potentially therapeutic agent.

25 In the context of the present invention the term to 'assess the effect' means to detect any changes in the severity or other characteristics of any one or more diseases in an individual. Such changes will be reflected in a change in sCD profile/fingerprint of an individual.

30 Preferably the agent is a potentially therapeutic agent. The term 'potentially therapeutic agent' means any agent that may cause a beneficial effect on an individual suffering from one or more diseases. Such beneficial effects may be for example reducing the clinical signs of the one or more diseases. It is an important feature of the

present invention though, that a change in level of any one sCD in isolation is not indicative of a change in severity of a disease.

- Generally, individual sCD levels will be elevated in a disease state as compared with a 'normal' non-diseased individual. Occasionally however, the level of an individual sCD will decrease in a disease state as compared with a normal non-diseased individual. However, according to the present invention, it is the changes of the profile of a number of sCDs (that is a fingerprint) during a disease which provides an accurate measure of the effect of one or more agents on a disease state in an individual.
- 10 One skilled in the art will appreciate that on occasion a selection of the complete repertoire of sCDs available for testing may be measured. The selection chosen will vary according to the disease state being tested.

- According to this aspect of the invention, sCDs suitable for generating a sCD fingerprint are as described herein.
- 15

- Therapeutic agents may be tested for their effect on any one or more disease states selected from the group consisting of the following: infections, autoimmune disease, neoplastic, vascular endocrinological, metabolic, inflammatory degenerative, psychiatric psychological, traumatic, drug/toxin-related, bacterial, fungal, protozoan and viral infections, non-neoplastic disorders; pain; diabetes, obesity; anorexia; bulimia; asthma; pregnancy; endocrine; vascular; metabolic; gastro-intestinal; iatrogenic; psychiatric; psychological; exercise-induced; diet-related; ME; degenerative; Parkinson's disease; thrombosis; atherosclerosis; acute heart failure; hypotension; hypertension; erectile dysfunction; urinary retention; metabolic bone diseases such as osteoporosis; angina pectoris; hepatitis; myocardial infarction; ulcers; allergies; rheumatoid arthritis; inflammatory bowel disease; irritable bowel syndrome benign prostatic hypertrophy; psychosis; psychiatric disorders; including anxiety; schizophrenia; manic depression; delirium; dementia; severe mental retardation and dyskinesias, such as Huntington's disease or Gilles de la Tourette's syndrome; and preferably tumours which can be benign or malignant cancers; breast cancer; myeloma; melanoma; bladder cancer; leukaemia; plasmacytoma and others, but most
- 20
- 25
- 30

preferably appendicitis; Bence Jones Proteinuria; Chronic Myeloid Leukaemia; Colorectal cancer; chronic renal failure; Crohn's Disease; Diabetic Nephropathy; Cardiac pathology; Infection; Liver damage; Lymphoma; macrocytic anaemia; Prostate Cancer; Oligoclonal Banding and PE/DVT

5

Suitable agents for assessment according to the method of the present invention may be naturally occurring or synthetic. Naturally occurring agents include proteins, peptides or nucleic acids. They may be agents known to be of therapeutic value or they may be of unknown therapeutic value.

10

In a further aspect, the present invention provides a method for sub-categorising a sCD fingerprint profile comprising the steps of identifying within one or more disease states one or more sub-group/s of sCDs wherein each sub-group of sCDs exhibits common characteristics distinguishing it from any other sub-group within that disease category.

15

As used herein the term a 'sCD sub-category' describes a sub-group of sCDs which show a defined fingerprint/profile (sub-fingerprint) of sCD levels within a larger fingerprint of one or more disease states wherein each sub-group of sCDs exhibits common characteristics distinguishing it from any other sub-group within those one or

20

more disease states.

In a final aspect the present invention provides a sCD reference database comprising pathological and/or normal sCD fingerprint patterns.

25

As herein described the term 'a database' refers to a collection of sCD fingerprints from normal 'non-diseased' and/or diseased individuals. Advantageously, the database is computer generated and/or stored.

Definitions

30

As herein described the term 'CD' refers to a different cell surface leukocyte molecule recognised by a given monoclonal or group of monoclonal antibodies which specifically 'cluster' to the antigen/molecule in question.

In the context of the present invention, the term 'secretory CD molecule (sCD)' refers to a secreted form of a cell surface leukocyte molecule recognised by a given monoclonal or group of monoclonal antibodies as herein described.

5

As used herein the term a 'sCD sub-category' describes a sub-group of sCDs which show a defined fingerprint/profile (sub-fingerprint) of sCD levels within a larger fingerprint of one or more disease states wherein each sub-group of sCDs exhibits common characteristics distinguishing it from any other sub-group within those one or more disease states.

10

As herein defined, the term 'secretory CD fingerprint (sCD)' describes the pattern or profile of levels of more than one secretory CD measured in one or more body fluids from one or more individuals. A sCD fingerprint as herein defined may be representative of one or more non-diseased individual or a one or more diseased individual/s. Preferably, a secretory CD fingerprint describes the level of five or more secretory CD molecules, more preferably it is 6, 7, 8, 9, 10 or more sCD molecules, more preferably still a secretory CD fingerprint describes the levels of 15 or more sCD molecules. Most preferably a secretory CD fingerprint describes the levels of 20 or more sCD molecules.

15

20

As defined herein, the term a 'disease state' refers to any impairment of the normal physiological functions affecting an organism or any disease condition, disorder or the presence of a particular microbial, viral, parasitic or other pathogenic agent known to one skilled in the art. Suitable disease states for analysis as described herein include but are not limited to: infectious, neoplastic, autoimmune, metabolic, degenerative, psychological, psychiatric, iatrogenic, inflammatory, drug or toxin related, vascular, traumatic and endocrine diseases. Advantageously, 'a disease state' as herein defined refers to any one or more disease selected from the group which includes but is not limited to: infections such as bacterial, fungal, protozoan and viral infections, non-neoplastic disorders; stroke; heart condition; atherosclerosis; pain; diabetes, obesity; anorexia; bulimia; asthma; Parkinson's disease; thrombosis; acute heart failure;

25

30

- hypotension; hypertension; urinary retention; metabolic bone diseases such as osteoporosis and osteopetrosis; angina pectoris; hepatitis; myocardial infarction; ulcers; asthma; allergies; rheumatoid arthritis; inflammatory bowel disease; irritable bowel syndrome benign prostatic hypertrophy; pancreatitis; chronic renal failure and
- 5 psychotic and neurological disorders, including anxiety, schizophrenia, manic depression, delirium, dementia, severe mental retardation and dyskinesias, such as Huntington's disease or Gilles de la Tourette's syndrome and others. Most preferably it refers to appendicitis; Bence Jones Proteinuria; Chronic Myeloid Leukemia; Colorectal cancer; chronic renal failure; Crohn's Disease; Diabetic Nephropathy;
- 10 Cardiac pathology; Infection; Liver damage; Lymphoma; macrocytic anaemia; Prostate Cancer; Oligoclonal Banding and Pulmonary Embolism/Deep Vein Thrombosis (eg DVT/PE). One skilled in the art will appreciate that this list is not intended to be exhaustive.
- 15 Examples of secretory cluster of differentiation molecules suitable for measurement to generate a sCD fingerprint for use in the methods of the present invention include but are not limited to CD14, CD25, CD31, CD44, CD50, CD54, CD62E, CD62L, CD86, CD95, CD106, CD116, CD124, CD138, CD141, CD40L, CD8, CD23, CD30, CD40. Those skilled in the art will be aware of other suitable sCD molecules for analyses
- 20 according to the methods of the present invention.

- As defined herein the term 'an antibody' includes within its scope for example IgG, IgM, IgA, IgD or IgE) or fragment (such as a FAb, F(Ab')₂, Fv, disulphide linked Fv, scFv, diabody) whether derived from any species naturally producing an antibody, or
- 25 created by recombinant DNA technology; whether isolated from serum, B-cells, hybridomas, transfectomas, yeast or bacteria).

Brief description of the figures

- 30 **Figure 1. Disease Groups. Multiples of upper limit of normal (ULN). All sCD's included**

The limits indicated by each point are:

- No shading $\leq 1 \times \text{ULN}$
- Lightly shaded $1-2 \times \text{ULN}$
- Darkly shaded $> 2 \times \text{ULN}$
- A white slash in the box indicates no data available.

5

Figure 2. All sCD's that appear not to discriminate from the normals (sCD's 21; 102; 117, 126; 130; 26; 44v5; 44v6; 62P).

10

Figure 3. Disease Groups. Mode of Response for Remaining 20 sCD's.
To simplify the data further the modal response for each disease group was plotted.

Figure 4. Disease Groups. Mode of Response for remaining sCD's.
Data has been ranked in order of increased expression.

Figure 5. Disease Groups. Mode of Response for remaining sCD's.

15

Figure 6. Disease Groups. Mode of Response for all sCD's.
As for Figure 3 (except all sCD's are included).
The limits indicated by each point are:

- No shading $\leq 1 \times \text{ULN}$
- Lightly shaded $1-2 \times \text{ULN}$
- Darkly shaded $> 2 \times \text{ULN}$

20

A white slash in the box indicates no data available.

Figure 7. Disease Groups. Mode of Response for all sCD's.

Figure 8. Disease Groups. Mode of Response for all sCD's.

25 Detailed description of the Invention
General techniques

Unless defined otherwise, all technical and scientific terms used herein have the same meaning as commonly understood by one of ordinary skill in the art (e.g., in cell culture, molecular genetics, nucleic acid chemistry, hybridisation techniques and biochemistry). Standard techniques are used for molecular, genetic and biochemical methods (see generally, Sambrook *et al.*, Molecular Cloning: A Laboratory Manual, 2d ed. (1989) Cold Spring Harbor Laboratory Press, Cold Spring Harbor, N.Y. and Ausubel *et al.*, Short Protocols in Molecular Biology (1999) 4th Ed, John Wiley & Sons, Inc. which are incorporated herein by reference) and chemical methods. In addition Harlow & Lane., A Laboratory Manual Cold Spring Harbor, N.Y, is referred to for standard Immunological Techniques.

Generation of a fingerprint of one or more disease states

In a first aspect, the present invention provides a secretory CD (sCD) fingerprint of one or more disease states.

15

Clinical signs and symptoms and various biochemical indicators of disease are used to identify individuals with one or more defined disease states. sCD levels are then measured for a number of sCDs present in one or more body fluid samples from each individual, preferably in a number of individuals using methods known to those skilled in the art and described herein, in order to generate a reference disease state or reference composite disease state fingerprint for those one or more given disease states.

20

(A) Diagnosis of disease states

25 i. Diagnostic Indicators Used.

Appendicitis

- Request for Amylase at admission A&E/MAU Subsequent Histopathological Diagnosis

30

Bence Jones Proteinuria

- Multiple myeloma in which the malignant plasma cells excrete only

light chains of one type (either 2 or 3); lytic bone lesions occur in about 60% of the cases, and light chains (Bence Jones protein) can be detected in the urine

- Positive finding

5

Chronic Myeloid Leukaemia

- Histopathological Diagnosis

Colorectal Carcinoma

10

- Histopathological Diagnosis

Chronic Renal Failure

- Prolonged elevation of serum creatinine.

15

Crohn's Disease

- Histopathological diagnosis

Diabetic Nephropathy

20

- Identified by abnormal urine Albumin/Creatinine ratio from subjects attending diabetic clinic.

Cardiac Pathology

- MI (as diagnosed by increased CK, symptoms and ECG changes).

25

Infection

CRP (C reactive protein) > 250 g/l (e.g. *Staphylococcus aureus* infections).

30

Liver Damage

- Clinical Details Alcoholic Liver Disease/Poisoning. Abnormal liver function tests.

Lymphoma

- Histopathological Diagnosis

5 Macrocytic Anaemia

- Diagnosed by haematological parameters. Hb <10 g/dL; MCV > 100 fL

Oligoclonal Banding

- 10 - small discrete bands in the gamma globulin region of the spinal fluid electrophoresis, indicating local central nervous system production of IgG; bands are frequently seen in patients with multiple sclerosis but can also be found in other diseases of the central nervous system including syphilis, sarcoidosis, and chronic infection or
- 15 inflammation.

VQ (T pulmonary angiogram)Pulmonary Embolism/Deep Vein Thrombosis

- 20 - ultrasound VQ or CT pulmonary angiogram scan (ventilation perfusion mismatch)

Prostate Carcinoma

- Histopathological diagnosis and elevated PSA.

25

In general, a combination of the patient's history, medical examination, general health and indicators provided from biochemical, histochemical, radiochemical and other types of tests and disease and/or clinical signs of disease will be used in the diagnosis of disease. For the avoidance of doubt, the term 'clinical signs and symptoms of

30 disease' means the same as 'clinical details' of disease.

(B) Samples of body fluids from disease states

For each sCD the following information is generally obtained a) the dynamic range of the assay b) the range of concentrations expected in health c) the range of concentrations expected in disease. From this information an approximate dilution factor for each assay may be obtained, allowing maximum use of subject samples. One skilled will appreciate thought that in some circumstances body fluid samples may not be diluted for testing.

Suitable body fluids for measuring sCD levels as herein defined include whole blood, serum, urine, tissue fluid, cerebrospinal fluid, lymphatic fluid, aspirate, bone marrow aspirate, mucus or other tissue or body fluid. One skilled in the art will appreciate that this list is not intended to be exhaustive. Preferably sCD levels are measured in serum which is prepared from whole blood using methods familiar to those skilled in the art. At least 1.5ml of sample is required for testing of all the sCDs. In addition, the present inventors have shown that haemolysis and lipaemia can interfere with some immunoassays used for detecting sCDs and therefore samples are used which exhibit minimal haemolysis and lipaemia.

Body fluid samples may be diluted in order to measure the sCD levels and the dilution factor for each sCD should be the same for the generation of the fingerprints for all disease states tested. One skilled in the art will appreciate that the dilution factor may be adjusted in order to focus on either high or low concentrations of sCDs. Advantageously, the dilution factor will be adjusted to focus on high concentrations of sCDs.

25 (C) Methods of measuring sCD levels.

Suitable methods for measuring levels of sCDs in body fluids include flow cytometry, in particular multiplexed particle flow cytometry, immunoassay and microarray technologies utilising antibody or ligand interactions. Advantageously, sCDs levels are measured using multiplexed particle flow cytometry and/or chip based monoclonal antibody technology. These methods will be familiar to those skilled in the art.

(i) Immunoassays

Immunoassays such as immunoblotting (detecting membrane-bound and soluble proteins), and enzyme linked immunoassays (ELISA) provide a sensitive and specific means of detecting target substances.

- 5 Although the various types of immunoassays are performed differently, they have one thing in common-they all involve antibodies. Used in an appropriate immunoassay system, specificity leads to sensitivity. As herein defined the term 'antibodies' includes antibody fragments, engineered immunoglobulin folds or scaffolds which have a binding affinity for soluble CD molecules.

10

One skilled in the art will appreciate that the 'immunoassay technique' may be adapted to use other molecules which selectively bind sCDs. Those skilled in the art will be aware of such molecules.

- 15 In a direct immunoassay, the antibody used as the primary reagent is advantageously given a fluorescent, enzymatic, or radio-active detection means. In indirect immunoassays, the secondary antibody-usually polyclonal antisera produced by a goat or a rabbit against human immunoglobulins-carries the detection means. When a secondary antibody is used, the initial immune reaction between the primary antibody and the target antigen is amplified, producing a more readily detectable signal.

20

- Western blots of electrophoretically separated proteins (immunoblots), on the other hand, are generally probed with antibodies labeled with an enzyme or a radioisotope such as ¹²⁵I. Chromogenic or chemiluminescent substrates can also be used. For example, enzymes such as HRP and AP catalyze chromagenic reactions, in which a colourless substrate is converted into a coloured compound, and also chemiluminescent reactions where light is emitted.

25

- Chromogenic substrate kits are commercially available and include but are not limited to for example alkaline phosphatase, horseradish peroxidase, and TMB peroxidase (TMB is tetramethylbenzidine, the substrate in this case). Boehringer Mannheim also has several enzyme substrates for immunoassays available. They include but are not

30

limited to for example ABTS (2,2'-azino-di-3-ethylbenzthiazoline sulfonate) and TMB (tetramethylbenzidine), which are used with HRP; and 4-nitrophenyl phosphate and 5-bromo-4-chloro-3-indolyl phosphate (BCIP) for immunoassays in which an alkaline phosphatase-conjugated antibody is used.

5

Chemiluminescent substrates are available from companies such as Pierce, which for example, produces the SuperSignal CL-HRP Substrate, an enhanced chemiluminescent substrate for horseradish peroxidase. This system detects specific proteins on immunoblots with a sensitivity that rivals radioactivity (reportedly to picogram levels).

10 When the chemiluminescent substrate is applied to membrane-bound proteins on an immunoblot, an instantaneous but long-lasting flash of light is produced.

Commercially available immunoassay kits for measuring sCD levels include those from Diaclone 1, Bd A. Fleming BP 1985 F-25020 Besancon Cedex-France which provides kits for the measurement of a number of CD molecules including CD 14, CD21, CD25, CD31, CD44, CD50, CD54, CD62E, CD62L, CD86, CD95, CD102, CD106, CD116, CD117, CD124, CD126, CD130, CD138, CD141, CD40L. Medsystems diagnostics GmbH, Rennweg 95b, A-1030 Vienna Austria, also provides kits which measure sCD levels.

20

(ii) Flow cytometry

Techniques for carrying out flow cytometry are familiar to those skilled in the art and are described in Flow Cytometry: A Practical Approach. Edited by MG Ormerod. IRL Press, Oxford. 1994. ISBN 0-19 963461-0. Practical Flow Cytometry. 3rd Edition.

25 Howard M Shapiro. Alan R Liss, Inc. ISBN 0-471-30376-3. Flow Cytometry. First Principles. Alice Longobardi Givan. Wiley-Liss, New York, 1992. ISBN 0-471-56095-2. Handbook of Flow Cytometry Methods. Edited by J Paul Robinson. Wiley-Liss, New York, 1993. ISBN 0-471-59634-5.

30 (iii) Multiplexed particle flow cytometry assay

Methods for simultaneously assaying different proteins in individual samples are commercially available. Some of those commercially available are detailed below:

The versatile laboratory multianalyte profiling (LabMAP™) system developed by Luminex Corp. of Austin, Texas, can be used for virtually any bioassay that is based on the specific binding of one molecule to another, for example a monoclonal antibody raised against a sCD and a CD molecule.

LabMAP assays for a sCD molecule can be based on the immunological detection, and/or may follow the gain or loss of fluorescence (e.g., when a mAb raised against a sCD binds to a sCD target). LabMAP assays employ three different fluorochromes: two to create color-coded microspheres, and the third for quantifying the reaction. Polystyrene microspheres are internally dyed with precise ratios of two spectrally distinct fluorochromes. This ratio confers a unique identifying "signature" or "spectral address" to each microsphere set.

Bioassays are conducted on the surfaces of the microspheres. Each capture probe (e.g., sCD-specific antibody or other affirmative detection reagent) is immobilized onto a color-coded set of microspheres using any of a variety of different surface chemistries. Luminex offers microspheres bearing Lumavidin™ (an avidin derivative), for immobilizing biotinylated molecules, or carboxyl groups, for covalently coupling protein. The binding of analyte to an immobilized probe is detected via a detection reagent labeled with the third fluorochrome. Luminex currently offers 100 different microsphere sets, each of which can be used for the simultaneous measurement of a different analyte. Thus, in theory, up to 100 different species can be simultaneously measured in a single tube or microplate well.

25

The Luminex microsphere product line is designed specifically to work with the instruments available from Luminex or their partners. The Luminex 100™ instrument uses microfluidics to align the microspheres in single file and employs two lasers, one for the detection of the fluorescent microsphere itself, and the other for the reporter reagent. The colour signals are captured by an optics system and translated into binding data via digital signal processing.

30

Instrumentation, reagents, and custom services for LabMAP technology users are also available from companies other than Luminex. For example, Bio-Rad Laboratories of Hercules, Calif., introduced the Bio-Plex™ Protein Array System. This system combines a fluorescent reader with the software, protocols, and supplies
5 needed for performing LabMAP-based assays in a 96-well microplate format. The primary benefits of the Bio-Plex Protein Array System, other than up to a 100-fold increase in data, include significantly reduced sample requirements.

The Cytometric Bead Array (CBA) system from BD Biosciences of San Diego is
10 flexible in that it accommodates multiple sizes and fluorescent intensities of particles. The system includes everything the researcher needs to implement this technology, including a cytometer setup kit with the requisite software, reagents and standards. The company's CBA assay kits employ their proprietary bead sets, which are internally
15 dyed with varying intensities of a proprietary fluorophore. These sets are distinguished via one fluorescence parameter and two size discriminators. However, the system is also capable of handling assays based on the use of other types of spectrally distinct microsphere sets. The CBA analysis software is an "add-in" for
20 Microsoft Excel®, and is compatible with contemporary data acquisition software such as CellQuest™. Researchers can employ a variety of preset configurations for generating standard dilution series and calibration curves, and data reports can be generated at each step in the process.

(iv) "Antibody Chip" array technology.

25 The array format has revolutionised biomedical experimentation and diagnostics, enabling ordered high-throughput analysis. During the past decade, classic solid phase substrates, such as microtitre plates, membrane filters and microscopic slides, have been turned into high-density, chip-like structures. The concept of the array library was central to this development which now extends from DNA to protein. Similar to
30 the gene chip arrays measuring RNA levels on a genome wide scale, recombinantly expressed genes can be robotically arrayed to produce protein chip arrays that can be used for the simultaneous assessment of protein levels on a proteome wide scale.

Additionally, protein specific antibodies can be arrayed to produce "antibody chip arrays" (Cahill D., 2001, *J. Imm. Meth.* 250, 81-91). The availability of such antibody chip arrays can be used to simultaneously analyse numerous interactions within a single sample. The "antibody chip" can be used to demonstrate antibody-protein interactions by incubating the chip with target proteins which have been labelled with a traceable marker (ProteinChip, CIPHERGEN Biosystems, Fremont, CA, USA; BIAcore chips, Biacore, Uppsala, Sweden) or by incubating the chip with protein molecules or fragments thereof and detecting association between antibody and protein molecule or fragment thereof using ELISA type assays (Caliper Technologies, Mountain View, CA, USA; Orchid Biocomputer Inc., Princeton, NJ, USA). In the context of the present invention, using antibody chip array technology is especially suitable as sCD molecules that can be found as complete or truncated/incomplete peptides will be recognised by a specific antibody. Chips may equally well be assayed with affirmative detection technologies.

15

Techniques for preparing antibody arrays are described below:

The antibodies may be covalently linked to a suitable membrane such as an Immobilon P membrane (PVDF; Millipore Corporation). Subsequent blocking with an excess of a protein solution such as a skim milk preparation is preferred. A blocking agent is designed to eliminate non-specific binding on the binding surface. Other suitable blocking agents are Irish moss extract or other source of carrageenan or gelatin. The antibodies are also adsorbed to a nitrocellulose film on a glass microscope slide (Schleicher and Schuell, NH, USA) and the unbound nitrocellulose is then blocked with skim milk. Antibodies are also adsorbed to Nylon membranes. To increase the accessibility of bound anti-CD antibodies to antigens on cells, the solid support used for the array is initially coated with a recombinant, truncated form of Protein G from *Streptococcus* which retains its affinity for the Fc portion of IgG lacks albumin and Fab binding sites, and membrane-binding regions (Goward et al., 1990). Antibodies are applied to this coat of Protein G and bind via their domains leaving the Fab domains free to interact with cells. The Fab domains are also further from the solid support providing greater accessibility of CD antigens on cell

membranes to antibodies.

The array of antibodies is also constructed on a membrane or a coverslip. In this case, the antibodies are covalently linked to the membrane as duplicate spots in a two
5 dimensional matrix. The spots are arranged in a matrix such as but not limited to a 15x 15 matrix.

The antibodies are advantageously monoclonal and are specific for the cluster of
differentiation (cluster designation) antigens (CD antigens). Details of CD antigens
10 are available at http://www.ncbi.nlm.nih.gov/prov/cd/index_molecule.htm. The spots are of microscopic size and are produced by the application of a drop (- 10 nanolitres) of antibody solution (e.g. 10/tg protein/ml) on designated portions of a membrane or glass surface such as a coverslip, first washed with a non-specific protein absorbent
15 such as 30% w/v skim milk (Dutch Jug, Bonlac Foods Ltd, Melbourne, Australia) and then rinsed. Other protein solutions and other brands of skim milk may also be employed. The antibodies may be covalently coupled to the solid support such as through amino groups of lysine residues, the carboxylate groups of aspartate or glutamic acid residues or the sulphydryl groups of cysteine residues. The array of
20 antibodies selectively binds cells from body fluids which express the respective antigens or may bind free antigens. A positive and/or negative control is included such as an antibody for surface molecules or soluble molecules known to be present in the sample. An example of one form of the assay device is shown in Fig. 3. The solid support is conveniently of similar size and shape to a microscope slide and may be constructed of glass or other polymeric material.

25

A wall around the microscope slide may be separately added or moulded with the slide and this facilitates retention of fluid material. The present invention extends to any other device capable of fulfilling the method of the present invention.

30 In the case of nitrocellulose based antibody arrays are preferably constructed using a Biodot Aspirate and Dispense System (Cartesian Technologies) where 5 nL dots are applied to a nitrocellulose film. on glass microscope slides (Schleicher and Schuell,

Cat. No. 10484182). Purified monoclonal antibodies (Beckman Coulter, Becton Dickinson or Biosource International) are used at concentrations recommended for flow cytometric analysis and are applied in the same buffers as supplied by the manufacturers. The nitrocellulose is then blocked by incubation with 5 % w/v skim milk (Dutch Jug) for 1.5 h at 37°C. These blocking conditions are chosen to minimize background binding.

The stability of the arrays is further enhanced by adding protein stabilizing agents to the antibodies (e.g. polyethylene glycol or stabilizer products commercially available from Surmodics, MN, USA).

The following description provides a preferred method for preparing the antibody arrays: The panel of antibodies is generally used to construct antibody arrays with the Cartesian Technologies PixSys™ 3200 Aspirate and Dispense System. The antibodies are chosen for use in a particular diagnosis or detection protocol. Each antibody is generally applied in the volume of from about 1 to about 10 nanoliters in a dot format at approximately from 0.5 to 1.5 min intervals to create an appropriate array on a nitrocellulose film generally laid on a solid support such as but not limited to a microscope slide. After dotting, the supports are assessed on a light box and the corners of the arrays marked gently using, for example, a lead pencil. The antibody arrays are then immersed in a blocking agent such as but not limited to skim milk, Irish moss extract or other source of carrageenan or gelatin. From about 2 % to about 15 % w/v skim milk powder in PBS at 4°C overnight or at 37°C for from about 60-120 minutes is particularly useful. After application of the blocking agent, the solid supports are washed gently with purified water and allowed to dry at room temperature for a period of time from about 60-120 minutes.

The solid supports are then stored in an airtight bag at 4°C in the dark

(D) A sCD fingerprint for one or more disease states according to the present invention
Advantageously, levels of sCDs are measured in diseased individuals and absolute values may be divided by the upper limit of normal (ULN) obtained from healthy

individuals. The data is collated and the resultant pattern of values obtained for each sCD for one or more given disease states, from one or more individuals forms the basis of a sCD fingerprint of one or more given disease states.

- 5 The statistical significance of the increases or decreases in sCD levels found in various disease states can be assessed using a number of methods.

The sCD fingerprint can advantageously be simplified by removing those sCDs whose levels do not generally change significantly during one or more given disease states.

- 10 Examples of such sCDs include but are not limited to CD21, CD102, CD117, CD 126, CD130, CD 26, CD44v5, CDv6, CD62P. For example see Figure 2.

- 15 To simplify the fingerprint further, instead of showing the data from each individual for any given sCD during a disease, a modal value for each sCD calculated from the group of individuals may be plotted. The rationale of this is demonstrated in Figures 3, 5, 6, 7, 8. This provides an easily readable and simple 'fingerprint' of a disease.

- One skilled in the art will appreciate that there are many methods suitable for the statistical analysis of the sCD level data measured as herein described. These include
20 but are not limited to cluster analysis and other statistical methods for the detection of patterns.

Uses of one or more 'fingerprints of disease' according to the present invention

- In a further aspect, the present invention provides a method for predicting the presence
25 of one or more disease states in an individual comprising the step of analysing the sCD fingerprint in that individual.

- In a further aspect still, the present invention provides a method for detecting the presence of one or more disease states in an individual comprising the step of
30 analysing the pattern/s of secretory CD levels of more than one secretory CD which is present in that individual.

In a further aspect still, the present invention provides a method for detecting the extent of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or more reference sCD fingerprint/s.

5

In yet a further aspect, the present invention provides a method for assessing the progression of a disease state in an individual comprising the step of comparing the sCD fingerprint of an individual at two or more periods during the occurrence of the disease.

10

In a further aspect still, the present invention provides a method for assessing the effect of one or more agent/s on one or more disease states in an individual comprising the step of comparing a sCD fingerprint of an individual at two or more time periods.

15

In a further aspect, the data generated by the present invention is used to compile a reference database, comprising pathological and/or normal sCD fingerprints, against which the expression sCD pattern of any individual will be compared.

20

In a further aspect still, the present invention provides the use of a sCD fingerprint to assess the effect of one or more agents in an individual.

An embodiment of the present invention is its use as a tool for assessing the affect different diet and exercise regimes may have on human or other mammals.

25

Additionally, the present invention may be used to construct sub-categories of sCD fingerprint profiles, suitable for common therapeutic treatment.

The invention will now be described by the following examples which are in no way limiting of the invention.

30

Example 1. Figure 1. Disease Groups. Multiples of ULN All sCD's included

Two values obtained (CD40L and CD30) for the individual, classified as normal, with a suspected drug overdose were omitted from the calculation of the upper limit of normal. The dilution factor for each sCD was fixed throughout the study. The results obtained are those of the diluted sample and have not been multiplied by the dilution factor. The absolute value of each data point was divided by the upper limit of normal (ULN) as defined above. Where the absolute value was greater than the dynamic range of the assay the result [9999] was recorded.

The limits indicated by each point are:

- Green $\leq 1 \times \text{ULN}$
- 10 • Blue $1 - 2 \times \text{ULN}$
- Red $> 2 \times \text{ULN}$
- A white block indicates no data available.

Example 2. Figure 2 Remove all sCD's that appear not to discriminate from the normals (sCD's 21; 102; 117; 126; 130; 26; 44v5; 44v6; 62P).

To simplify the diagram the above sCD plots were removed.

The data suggests (Figure 1.) that the concentration of some of these sCD may actually be lowered in disease. As we initially worked on the premise that there would be over-expression of these molecules in disease, samples have been diluted optimally to focus on high, rather than low concentrations.

Example 3. Figure 3. Disease Groups. Mode of Response for Remaining 20 sCD's

To simplify the data further the modal response for each disease group has been plotted. As the lymphoma and "oligoclonal-banding positive" group contain only a single subject, they have been omitted. Where there is no clear mode, both responses have been shown.

Example 4. Figure 4. Disease Groups. Mode of Response for remaining sCD's

Data has been ranked in order of increased expression.

Example 5. Figure 5. Disease Groups. Mode of Response for remaining sCD's

As for Figure 4, except responses have been classified as "normal" and "abnormal".
Values >1 ULN have been classified as abnormal.

5

Both Figures 4 and 5 suggest that each disease state exhibit a unique pattern of elevated sCD expression.

Example 6. Figure 6. Disease Groups. Mode of Response for all sCD's

10

As for Figure 3 (except all sCD's are included).

The limits indicated by each point are:

- Green $\leq 1 \times \text{MoM}$
- Blue $1-2 \times \text{MoM}$
- 15 • Red $> 2 \times \text{MoM}$
- A white block indicates no data available.

Example 7. Figure 7. Disease Groups. Mode of Response for all sCD's

As for Figure 4.

Example 8. Figure 8. Disease Groups. Mode of Response for all sCD's

20

As for Figure 5. (except that values > 2 MoM have been classified as abnormal).

Comparison of figures 5 and 8 illustrate the importance of determining the cut-off threshold values in order to obtain a defined pattern.

25

All publications mentioned in the above specification are herein incorporated by reference. Various modifications and variations of the described methods and system of the present invention will be apparent to those skilled in the art without departing from the scope and spirit of the present invention. Although the present invention has

been described in connection with specific preferred embodiments, it should be understood that the invention as claimed should not be unduly limited to such specific embodiments. Indeed, various modifications of the described modes for carrying out the invention which are obvious to those skilled in biochemistry, molecular biology
5 and biotechnology or related fields are intended to be within the scope of the following claims.

CLAIMS

1. A secretory CD (sCD) fingerprint of one or more disease states.
- 5 2. A method of generating a secretory CD (sCD) fingerprint of one or more disease state/s comprising the step of measuring the levels of more than one secretory CDs from one or more individuals and collating the data.
- 10 3. A sCD fingerprint according to claim 1 or a method according to claim 2 wherein the disease state is any one or more selected from the group consisting of: infectious, neoplastic, autoimmune, metabolic, degenerative, psychological, psychiatric, iatrogenic, inflammatory, drug or toxin related, vascular, traumatic and endocrine diseases.
- 15 4. A sCD fingerprint or a method according any preceding claim wherein the disease state is any one or more selected from the group consisting of the following: infection, Bence Jones Proteinuria, Chronic Myeloid Leukemia, Colorectal cancer, chronic renal failure, Crohn's Disease, Diabetic Nephropathy, Cardiac pathology, Infection, Liver damage, Lymphoma, macrocytic anaemia, Prostate Cancer, Oligoclonal Banding and Pulmonary
20 Exubolism/Deep Vein Thrombosis and appendicitis.
- 25 5. A sCD fingerprint according to claim 1 or claim 3 or claim 4 or a method according to claim 2, claim 3 or claim 4 wherein the sCDs referred to comprise two or more selected from the group consisting: CD14, CD25, CD31, CD44, CD50, CD54, CD62E, CD62L, CD86, CD95, CD106, CD116, CD124, CD138, CD141, CD40L, CD8, CD23, CD30, CD40 and their homologues present in other mammalian or non-mammalian species.
- 30 6. A method according to any of claims 2 to 5 wherein the sCD levels are measured in samples of one or more body fluids from an individual.

35

7. A method according to claim 6 wherein the body fluid is serum.
8. A method according to any of claims 2 to 7 wherein sCD levels are measured using one or more methods selected from the group consisting of:
5 immunoassay and flow cytometry.
9. A method according to claim 8 wherein sCD levels are measured using multiplexed particle flow cytometry and/or chip based monoclonal antibody technologies.
10
10. A method for predicting the presence of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or more reference sCD fingerprint/s.
- 15 11. A method for detecting the presence of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated from that individual with one or more reference sCD fingerprint/s.
12. A method for detecting the extent of one or more disease states in an individual comprising the step of comparing one or more sCD fingerprint/s generated
20 from that individual with one or more reference sCD fingerprint/s.
13. A method for assessing the progression of a disease state in an individual comprising the step of comparing the sCD fingerprint of an individual at two or
25 more periods during the life-span of the disease.
14. A method for assessing the effect of one or more agent/s on one or more disease states in an individual comprising the step of comparing a sCD fingerprint of an individual at two or more different time periods.
30
15. The use of a sCD fingerprint to assess the effect of one or more agent/s on an individual.

16. A method for sub-categorising a sCD fingerprint profile comprising the steps of identifying within one disease category one or more group/s of sCDs wherein each group of sCDs exhibits common characteristics distinguishing it from any other group within that disease category.
- 5
17. A sCD database comprising pathological and/or normal sCD fingerprint patterns.
- 10

ABSTRACT

The present invention relates to the use of cluster of differentiation (CD) molecules in detecting the presence and progression of one or more disease states in an individual. In particular it relates to the use of profiles of secretory CD (sCD) molecules in
5 detecting and assessing the progression of one or more disease states in an individual. Further uses of sCD profiles according to the present invention are also described.

0030819 07-Mar-02 04:44

	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD	CD
	14	25	31	44	50	54	62	62	86	95	106	116	124	138	141	40	8	23	30	40	
Appendicitis	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Bence Jones Proteinuria	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Chronic Myeloid Leukaemia	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Colorectal Cancer	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Chronic Renal Failure	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Crohn's Disease	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Diabetic Nephropathy	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
High CK	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Infection	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Liver damage	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Low CK	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Lymphoma	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	
Macrocytic Anaemia	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0			

07-03-02

04:18pm

From-D YOUNG & CO

+02380 719800

T-511

P.044/161

F-513

Figure 3.

Disease	CD14	CD23	CD31	CD44	CD50	CD64	CD62E	CD62L	CD68	CD25	CD108	CD115	CD124	CD138	CD141	CD40L	CD8	CD23	CD31	CD47
Appendicitis	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
DJP	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
CML	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Cholangiocarcinoma	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
CRP	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Crohn's	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
DM	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
High CK	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Infection	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Liver Damage	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Low CK	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
MA	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
PCPVT	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Prostate	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Normal	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
CCI	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
CC2	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1

Figure 4.

Disease	CD25	CD30	CD134	CD140L	CD143	CD54	CD14	CD138	CD62E	CD44	CD126	CD114	CD22	CD40	CD8	CD56	CD300	CD9	CD62L	CD68
CRP																				
Liver Damage																				
MA																				
Infection																				
Pept V																				
High CK																				
Low CK																				
Appendicitis																				
BIP																				
Colorectal																				
CML																				
DM																				
Cytoma																				
Prostate																				
Normal																				
CD7																				
CD2																				

07-03-02

04:20pm

From-D YOUNG & CO

+02380 719800

T-511 P.046/161 F-518

Figure 5.

Disorder	CD25	CD141	ED108	CD30	CD54	CD40L	CD14	CD62E	CD44	CD134	CD8	CD138	CD85	CD50	CD11b	CD22	CD31	CD62L	CD40	CD98
CRF													1 1					1 1		1 1
Liver Damage													1 1	1 1	1 1			1 1	1 1	1 1
PEIOVE				1 1							1 1	1 1			1 1	1 1			1 1	1 1
CML	1 1	1 1											1 1	1 1	1 1	1 1			1 1	1 1
MIA							1 1						1 1	1 1	1 1	1 1		1 1	1 1	1 1
Infection			1 1										1 1	1 1	1 1	1 1			1 1	1 1
Appendicitis				1 1									1 1	1 1	1 1	1 1			1 1	1 1
BAP													1 1	1 1	1 1	1 1			1 1	1 1
Colorectal													1 1	1 1	1 1	1 1			1 1	1 1
Low CK													1 1	1 1	1 1	1 1			1 1	1 1
High CK	1 1	1 1											1 1	1 1	1 1	1 1			1 1	1 1
Cervix		1 1		1 1									1 1	1 1	1 1	1 1			1 1	1 1
Prostate				1 1									1 1	1 1	1 1	1 1			1 1	1 1
DN	1 1												1 1	1 1	1 1	1 1			1 1	1 1
Normal	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1
CD3	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1	1 1
CCR																				

07-09-02

04:20pm

From-D YOUNG & CO

+02380 719800

T-511 P.047/161 F-513

Figure 6.

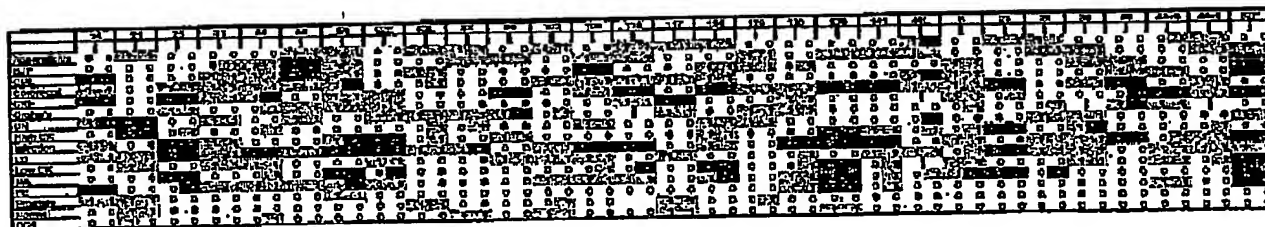


Figure 8.

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32	33	34	35	36	37	38	39	40	41	42	43	44	45	46	47	48	49	50	51	52	53	54	55	56	57	58	59	60	61	62	63	64	65	66	67	68	69	70	71	72	73	74	75	76	77	78	79	80	81	82	83	84	85	86	87	88	89	90	91	92	93	94	95	96	97	98	99	100	101	102	103	104	105	106	107	108	109	110	111	112	113	114	115	116	117	118	119	120	121	122	123	124	125	126	127	128	129	130	131	132	133	134	135	136	137	138	139	140	141	142	143	144	145	146	147	148	149	150	151	152	153	154	155	156	157	158	159	160	161	162	163	164	165	166	167	168	169	170	171	172	173	174	175	176	177	178	179	180	181	182	183	184	185	186	187	188	189	190	191	192	193	194	195	196	197	198	199	200	201	202	203	204	205	206	207	208	209	210	211	212	213	214	215	216	217	218	219	220	221	222	223	224	225	226	227	228	229	230	231	232	233	234	235	236	237	238	239	240	241	242	243	244	245	246	247	248	249	250	251	252	253	254	255	256	257	258	259	260	261	262	263	264	265	266	267	268	269	270	271	272	273	274	275	276	277	278	279	280	281	282	283	284	285	286	287	288	289	290	291	292	293	294	295	296	297	298	299	300	301	302	303	304	305	306	307	308	309	310	311	312	313	314	315	316	317	318	319	320	321	322	323	324	325	326	327	328	329	330	331	332	333	334	335	336	337	338	339	340	341	342	343	344	345	346	347	348	349	350	351	352	353	354	355	356	357	358	359	360	361	362	363	364	365	366	367	368	369	370	371	372	373	374	375	376	377	378	379	380	381	382	383	384	385	386	387	388	389	390	391	392	393	394	395	396	397	398	399	400	401	402	403	404	405	406	407	408	409	410	411	412	413	414	415	416	417	418	419	420	421	422	423	424	425	426	427	428	429	430	431	432	433	434	435	436	437	438	439	440	441	442	443	444	445	446	447	448	449	450	451	452	453	454	455	456	457	458	459	460	461	462	463	464	465	466	467	468	469	470	471	472	473	474	475	476	477	478	479	480	481	482	483	484	485	486	487	488	489	490	491	492	493	494	495	496	497	498	499	500	501	502	503	504	505	506	507	508	509	510	511	512	513	514	515	516	517	518	519	520	521	522	523	524	525	526	527	528	529	530	531	532	533	534	535	536	537	538	539	540	541	542	543	544	545	546	547	548	549	550	551	552	553	554	555	556	557	558	559	560	561	562	563	564	565	566	567	568	569	570	571	572	573	574	575	576	577	578	579	580	581	582	583	584	585	586	587	588	589	590	591	592	593	594	595	596	597	598	599	600	601	602	603	604	605	606	607	608	609	610	611	612	613	614	615	616	617	618	619	620	621	622	623	624	625	626	627	628	629	630	631	632	633	634	635	636	637	638	639	640	641	642	643	644	645	646	647	648	649	650	651	652	653	654	655	656	657	658	659	660	661	662	663	664	665	666	667	668	669	670	671	672	673	674	675	676	677	678	679	680	681	682	683	684	685	686	687	688	689	690	691	692	693	694	695	696	697	698	699	700	701	702	703	704	705	706	707	708	709	710	711	712	713	714	715	716	717	718	719	720	721	722	723	724	725	726	727	728	729	730	731	732	733	734	735	736	737	738	739	740	741	742	743	744	745	746	747	748	749	750	751	752	753	754	755	756	757	758	759	760	761	762	763	764	765	766	767	768	769	770	771	772	773	774	775	776	777	778	779	780	781	782	783	784	785	786	787	788	789	790	791	792	793	794	795	796	797	798	799	800	801	802	803	804	805	806	807	808	809	810	811	812	813	814	815	816	817	818	819	820	821	822	823	824	825	826	827	828	829	830	831	832	833	834	835	836	837	838	839	840	841	842	843	844	845	846	847	848	849	850	851	852	853	854	855	856	857	858	859	860	861	862	863	864	865	866	867	868	869	870	871	872	873	874	875	876	877	878	879	880	881	882	883	884	885	886	887	888	889	890	891	892	893	894	895	896	897	898	899	900	901	902	903	904	905	906	907	908	909	910	911	912	913	914	915	916	917	918	919	920	921	922	923	924	925	926	927	928	929	930	931	932	933	934	935	936	937	938	939	940	941	942	943	944	945	946	947	948	949	950	951	952	953	954	955	956	957	958	959	960	961	962	963	964	965	966	967	968	969	970	971	972	973	974	975	976	977	978	979	980	981	982	983	984	985	986	987	988	989	990	991	992	993	994	995	996	997	998	999	1000
1A																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																								

THE PATENT OFFICE

3 APR 2003

Received in Patents
International Unit

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

☐ BLACK BORDERS

☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES

☒ FADED TEXT OR DRAWING

☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING

☐ SKEWED/SLANTED IMAGES

☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS

☐ GRAY SCALE DOCUMENTS

☐ LINES OR MARKS ON ORIGINAL DOCUMENT

☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY

☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.